

Validation of novel *CNBP* repeat features in Myotonic Dystrophy Type 2 through comparative long-read sequencing with Oxford Nanopore and PacBio

^{1,2}Suárez-López FJ, ¹Carlomagno M., ¹Maestri S., ¹Esposito A., ³Obadović V., ³Ciabini D., ⁴D'apice M.R., ⁴Novelli G., ^{1,2}Delledonne M., ³Botta A., ^{1,2}Rossato M.

¹Department of Biotechnology, University of Verona, Verona, Italy
²Genartis s.r.l., Via Albere, Verona, Italy
³Department of Biomedicine and Prevention, University of Rome Tor Vergata, Rome, Italy
⁴Laboratory of Medical Genetics, Tor Vergata Hospital, Rome, Italy



To put the correct code



UNIVERSITÀ di VERONA




Introduction


Myotonic Dystrophy Type 2 (DM2) is caused by massive and highly mosaic CCTG repeat expansions in the *CNBP* gene. While targeted long-read sequencing has enabled the sequencing of *CNBP* repeats and revealed a novel TCTG interruption pattern (Alfano M. et al, eLife 2022), the computational characterization of this region remains challenging. Currently, there is a lack of clear consensus and benchmarking of bioinformatic pipelines to accurately identify the architecture of these expansions. **In this work, our goal was to establish a standardized computational framework through systematic algorithmic benchmarking and cross-platform validation, in order to accurately resolve these expansions**

Methods


1. Targeted sequencing using LRS

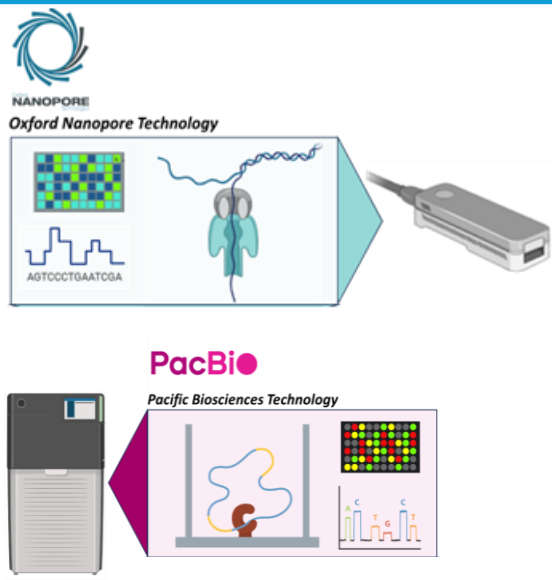


2. Basecaller benchmarking (ONT)

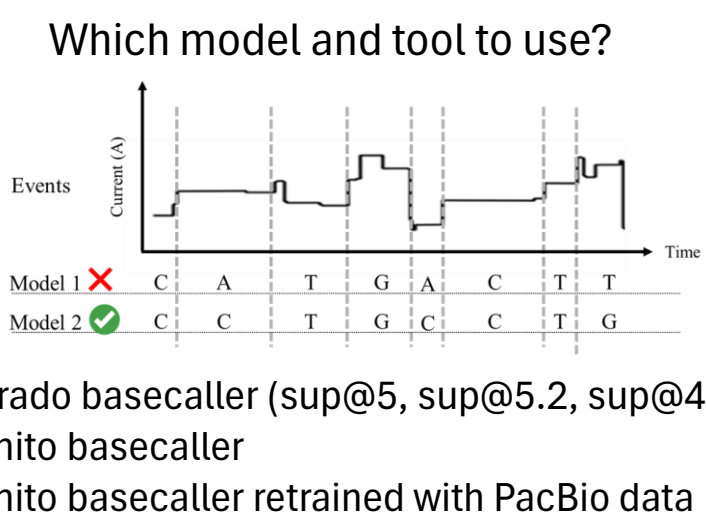

3. Downstream analysis

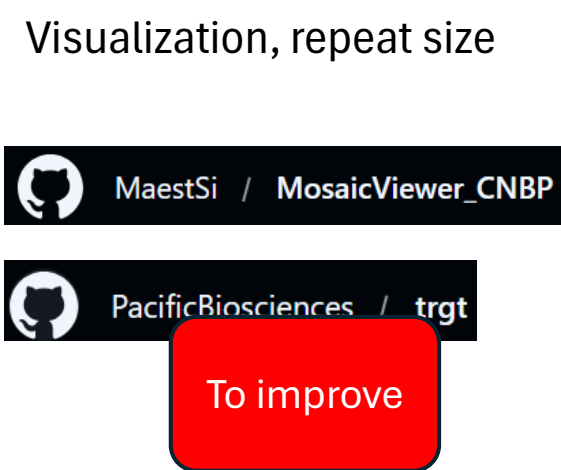


4. Training Bonito basecaller to improve the accuracy

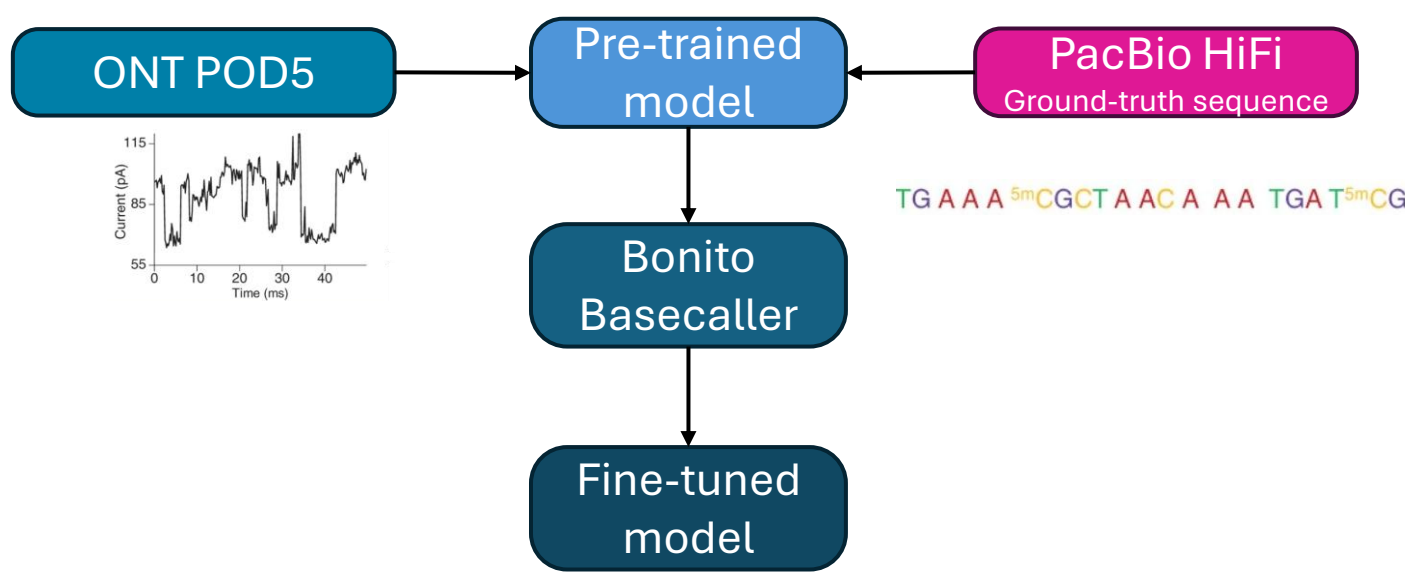


We sequenced the largest DM2 cohort to date, using a Cas9-mediated approach combined with Oxford Nanopore (ONT) and PacBio platforms. To assess the impact of signal decoding on repeat resolution, raw POD5 files from ONT data were basecalled using multiple versions of the Dorado super-accurate model. Complete reads spanning the entire expansion were extracted for downstream benchmarking. **We compared two distinct tools for expansion sizing, motif identification, and mosaicism visualization: TRGT and MosaicViewer, developed for repeat expansion analysis with PacBio and ONT, respectively.**

In order to improve the accuracy of the models when basecalling highly repetitive motifs, we fine-tuned a pre-trained model using Bonito Basecaller, using the ONT POD5 data and, as their ground-truth reference, the HiFi sequences obtained with PacBio


Results

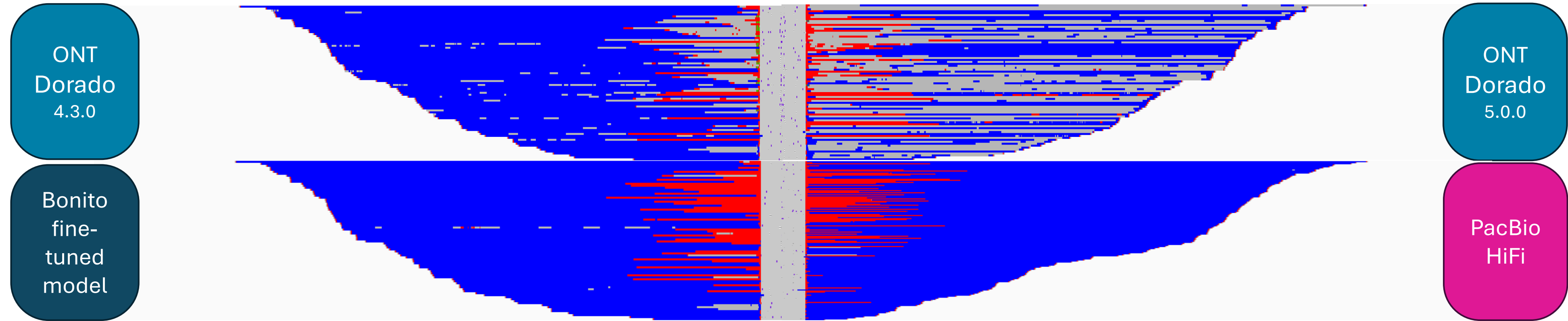


Our benchmarking revealed that the choice of the basecalling algorithm version significantly impacts the accuracy of the repeat sequence identification based on ONT data. In particular, specific Dorado versions showed higher concordance with PacBio HiFi data, proving that basecaller version has a great impact on the resolution of repeat motif

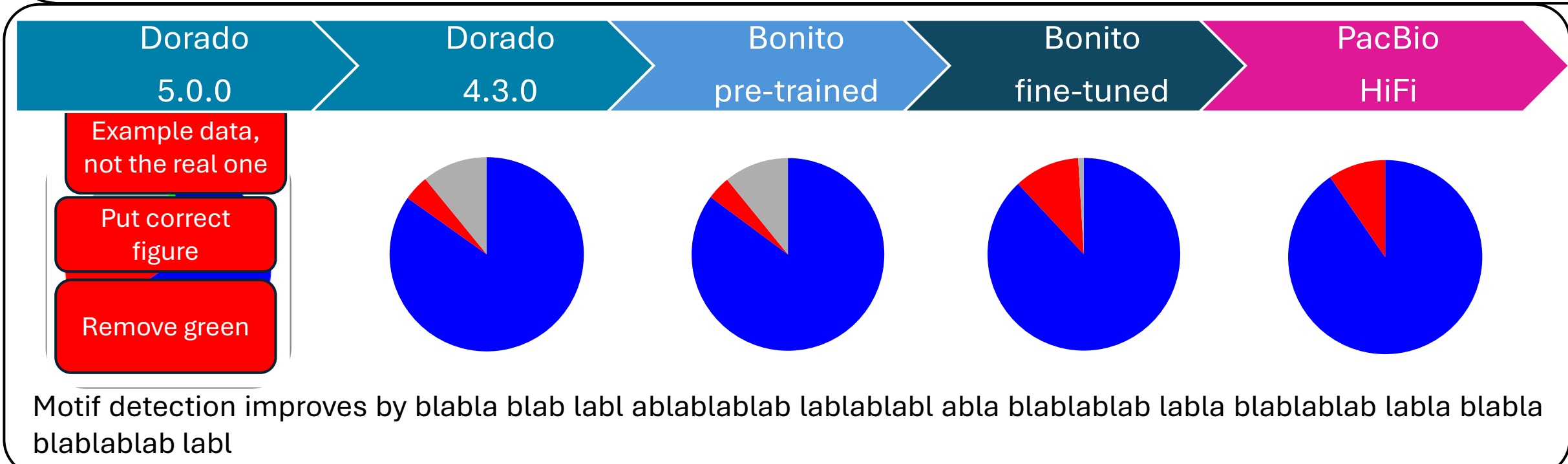

Motif detection improves by blabla blab labl ablablablab lablablabl abla blablablab labla blablablab labla blabla blablablab labl

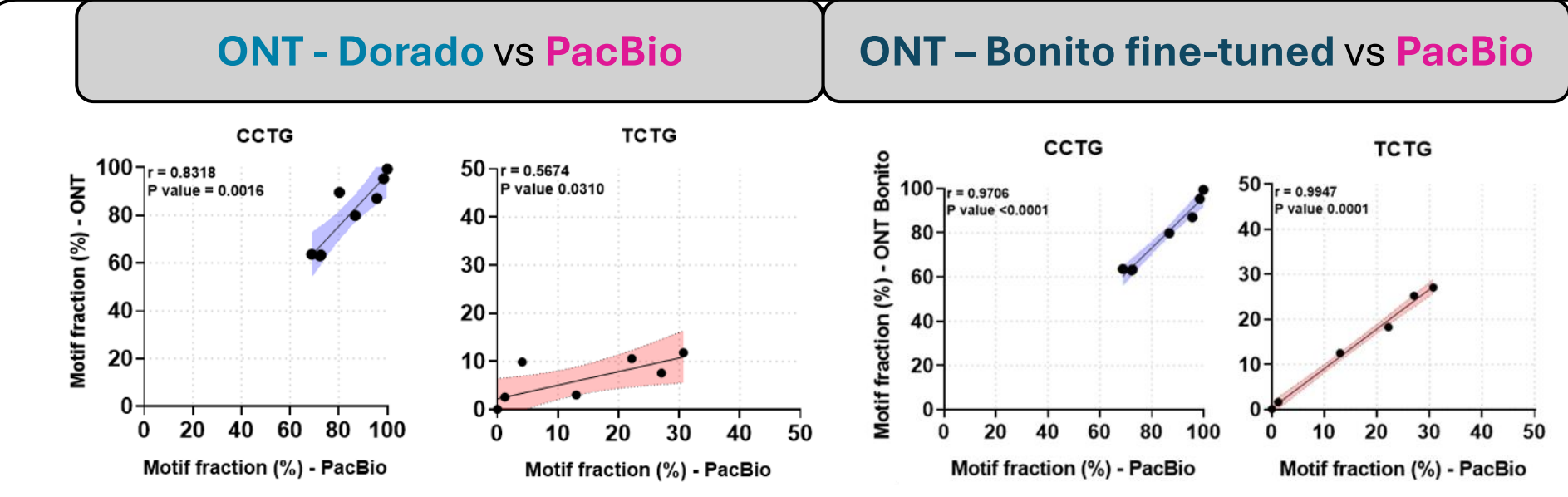

Correlations between the fraction of repeat units across the total number of repeats show an improvement after fine-tuning Bonito

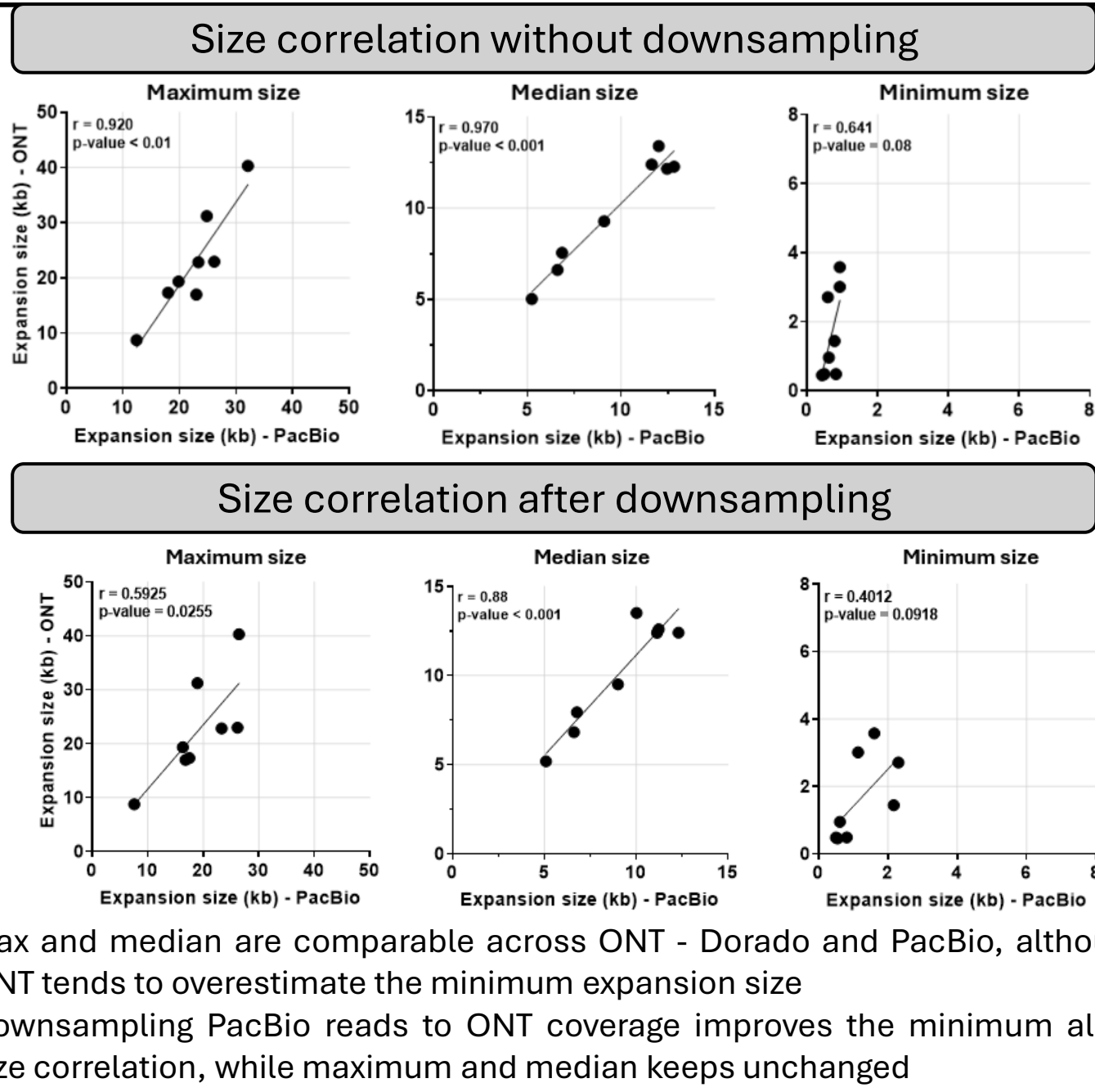

Max and median are comparable across ONT - Dorado and PacBio, although ONT tends to overestimate the minimum expansion size
Downsampling PacBio reads to ONT coverage improves the minimum allele size correlation, while maximum and median keeps unchanged

Conclusions

1. The basecaller impacts on the accuracy...
2. Dorado is comparable with PacBio for every aspect (size, mosaicism, interruptions...)
3. Novel TCTG motif
4. TRGT is comparable with MosaicViewer.....
5. Bonito ...